


SYRINGE PARTS	CLASSIFICATION
	<p>The disposable syringe with needle pertains to IIa class for medical devices, conforming to Rule 6 Enclosure IX directive 93/42 CEE of the Council June 14, 1993.</p> <p style="text-align: center;">CE Mark: 0123</p>
	MANUFACTURER
	<p>Pentaferte s.r.l. S.P. 262 Località Nocella 64012 Campi (TE) Phone: +39.0861.560201 Fax: +39.0861.560202</p>
	CND CODE
	A020102020101
GMDN CODE	
47017	

WARNING
<ul style="list-style-type: none"> • The CRN® CRONO® reservoirs have been designed to be used exclusively with CRONO® model pumps. • The use of CRN® CRONO® reservoirs with other pumps, not of the Crono series, can cause damage to the pump and injury to the patient. • CANÈ S.p.A. declines all responsibility if: <ul style="list-style-type: none"> - CRN® CRONO® reservoirs are used with other devices - CRONO® model pumps are used with reservoirs other than original CRN® CRONO® reservoirs

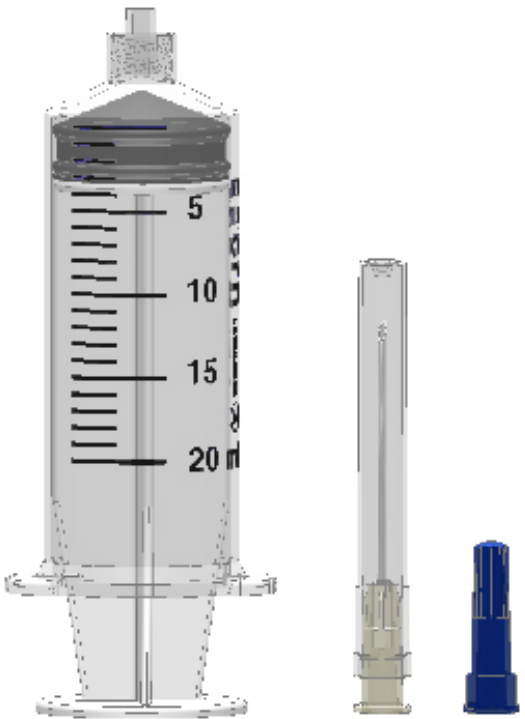

CRN® CRONO® SYRINGE 20 ml SPECIFICATIONS

Rev. 10 of 01.12.2014

COMPONENT (Drawing's code)	FEATURES (Standard's reference)	MATERIAL (Data sheet)
20ml rod (CRN/01b/B)	Flange cylinder/flange distance: min.20mm (ISO 7886-2) The head is threaded to unscrew the rubber piston	Matt Polypropylene medical degree HF136MO
20ml cylinder (CRN/03/D)	Utilizable capacity: min. 22 ml Luer Lock cone (ISO 594/1 and 594/2)	Transparent polypropylene medical degree RF825MO
	Scale Graduation space: 1ml (ISO 7886-1) Scale numeration: every 5ml (ISO 7886-1) Capacity output tolerance beyond 5ml: +/-4% (ISO 8537)	Black ink
	Lubrication: maximum quantity 0,25 mg/cm2 of internal surface (ISO 7886-1)	Silicone oil + isooctane diluent
20ml piston (CRN/05/C)	Black colour; two hold rings Lubricant: silicone viscosity 1000 Internal screw threading with head/rod	Latex and 2MBT Free rubber
Needle 19 G x 11/2" on side	Barrel: cone Luer 6% (ISO 594/1) Standard coloring ISO	Transparent polypropylene medical degree + coloured master medical degree
	Needle cover	Transparent polypropylene medical degree
	Cannula: Triple sharpening with primary angle 11°+/-2° (ISO 7864) External diameter and internal diameter complying to ISO 9626 standard Length complying to ISO 7864 standard Cleaning: electrocleaning + ultrasounds	Stainless steel AISI 304 Epoxy resin Silicone Polydymetilsiloxane (complying to ISO 7864 and FDA standard)
Luer-lock cap (CRN/25)	Blue cap: cone Luer 6% (ISO 594/1) to inserting on Luer Lock cone	Polypropylene medical degree + coloured master medical degree
Assembled syringe (CRN/04/C)	Hold while compressing: min. 300 kPa (ISO 7886-1) Dead space: max. 0,10ml (ISO 7886-1)	-----

COMPONENT (Drawing's code)	FEATURES (Standard's reference)	MATERIAL (Data sheet)
Packing	Paper: 80 g/m2 Stamping: Product description, Lot Nr., Expiring date	Medical degree paper Film
	Film: 90 µm thickness Type of opening/blister: PEEL	PP/PE
Box	Number of pieces per dinking die: 80 External dim.: (mm.) 382 x 232 x 183 Stamping: Product description, Lot Nr., Expiring date	Corrugated cardboard with white kraft external cover
Master	Number of dinking die per master: 4 (total 320 pcs) External dim.: (mm.) 394 x 394 x 482 Stamping: Product description, Lot Nr., Expiring date, Sterilization's cycle Nr.	Two waves microtriple corrugated cardboard with white kraft external cover
Sterilization	Mixture: Ethylene Oxide 10%+CO2 90% Autoclaves capacity: 24 m3 and 40 m3 Cycle validation notified by external SINAL accredited Laboratory, according to EN 550 standard "Medical devices' sterilization - Routine's validation and control of ethylene oxyde sterilization"	

EXPIRY DATE (SHELF LIFE)	STORAGE TEMPERATURE
5 years from the manufacture date	-10°C to 40°C

CANÈ S.p.A. Tecnologie Medicali - Via Cuorgnè 42/a - Rivoli (TO) Tel. 011.9574872 - Fax 011.9598880
Internet: www.canespa.it E-mail: mailbox@canespa.it